

# H2O<sub>A</sub> Hypoxia workstation

Exclusively at:  
**HYPOXYGEN**  
v. Define your environment

## Lift-off Lid Facility

The transparent top provides excellent visibility of the working and incubation areas. The top lifts off for deep cleaning, full sample ingress/egress and the introduction of laboratory instrumentation.

## Optional Internal Socket

Some users need to operate electrical instrumentation inside the workstation. A fused mains power supply and country-specific electrical socket is available as an option.

## Tissue Flask Rack

Bespoke racks to accommodate popular sizes of cellware.

## Refillable Chemical Sachets

Catalyst for oxygen reduction and patented chemical compounds for the removal of hydrogen sulphide and volatile fatty acids are introduced when anoxic conditions are required.

## Multi-functional Porthole System

The unique, patented, oval portholes provide both sample transfer and operator entry and exit.

Each porthole can be used to transfer up to 6 x 50ml flasks, or 2 x 275ml flasks or 1 x 600ml flask (or objects of similar size) at the same time as an operator's arms are inserted or withdrawn from the workstation.


## Automatic Humidity Control System

The fully automated humidity control system does not require any user maintenance.

## Capacity

A space measuring 655mm wide x 255mm deep x 500mm high (25.75" x 10" x 19.5") can accommodate samples and equipment whilst still providing a generous working area.

Adjustable height shelving and bespoke racking are available as options.

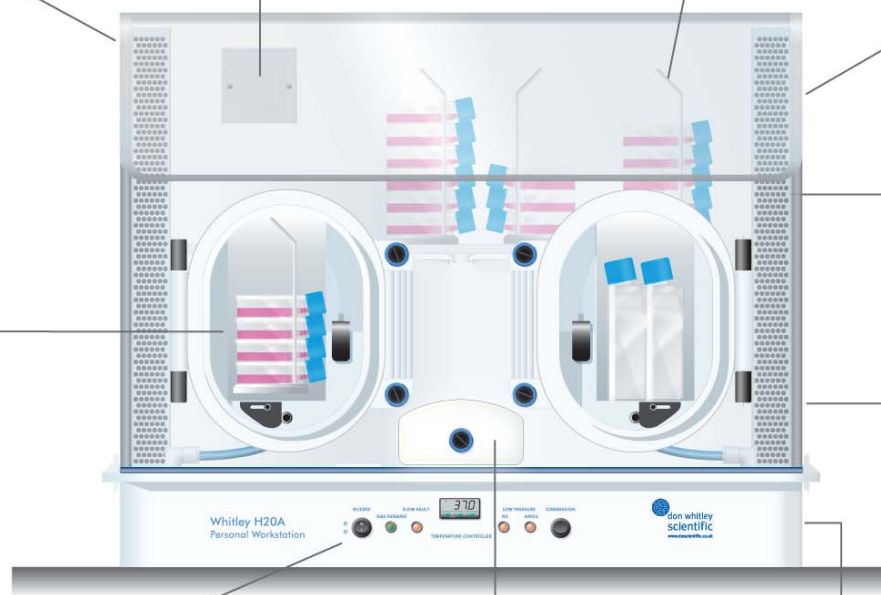

## 10 Year Shell Life

Manufactured from 10mm acrylic sheet and annealed during fabrication for maximum strength and durability.

## Easy Access to all Control Systems

All components are easily accessed from the front of the workstation.

## Single Sample Entry System

Individual 50ml flasks and 90mm Petri dishes can be introduced quickly and conveniently into the workstation.

## Integrated Handles

Two handles each side make workstation movement and lifting easy.

## Adjustable Height Trolley

An adjustable height trolley is available as an option.

877-HYPOXYGEN (497-6994)  
www.HYPOXYGEN.us