

# VACUBOY



Multifunctional Hand Operator for Vacuum Systems

**HYP**<sub>2</sub>**OXYGEN**  
v. Define your environment

# VACUBOY

## The VACUBOY fits all...

The **VACUBOY** hand operator is designed to accommodate a variety of aspiration adapters – either stainless steel, Pasteur/serological pipettes or disposable tips. Each adapter is designed to perfectly match the physical dimensions of typical laboratory glass and plastic containers.

The **VACUBOY** is a very **convenient tool** to control aspiration and disposal of fluids from slides, cuvettes, tubes, micro plates, Petri-dishes, T-flasks, roller bottles, bioreactors and other typical containers used in laboratories.

Furthermore it fits any existing vacuum source in a laboratory, therefore it is an **ideal upgrade** to the infrastructure already in place.

### Benefits

- Very easy and time saving operation
- Ergonomic design takes the strain out of tedious aspiration tasks

### Customers

- Pharmaceutical / Biotech industry
- Cell culture facilities
- Laboratories working with bio-hazardous fluids
- Laboratories working with pathogenic material

### Features


Single and Multi-channel aspiration adapters:

- Single channel
- 4-channel
- 8-channel

For different tips:

- Pasteur/serological pipettes
- Stainless steel
- Disposable tips

All parts in contact with fluids are autoclavable

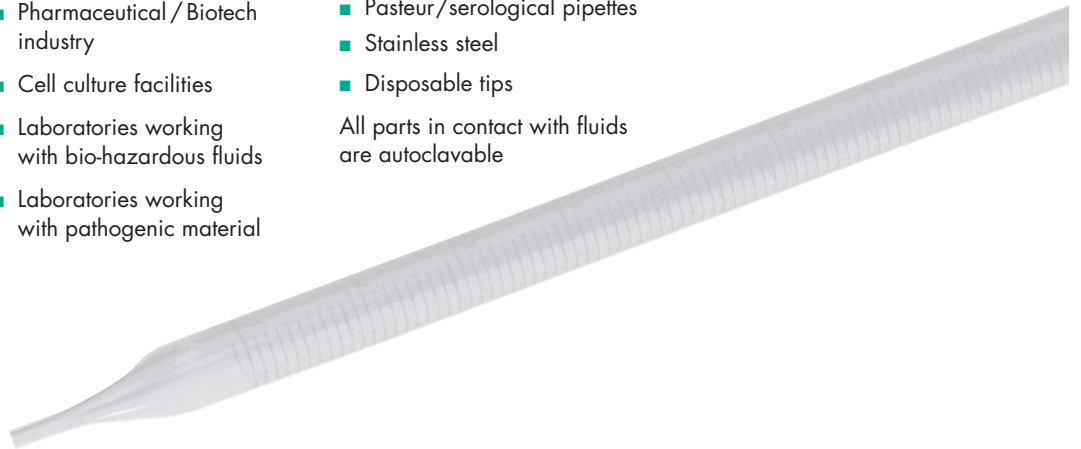

## ...Applications



T-Flasks



Roller Bottles



Microplates




Erlenmeyer Flasks



6-well Plates


...Vacuum Sources

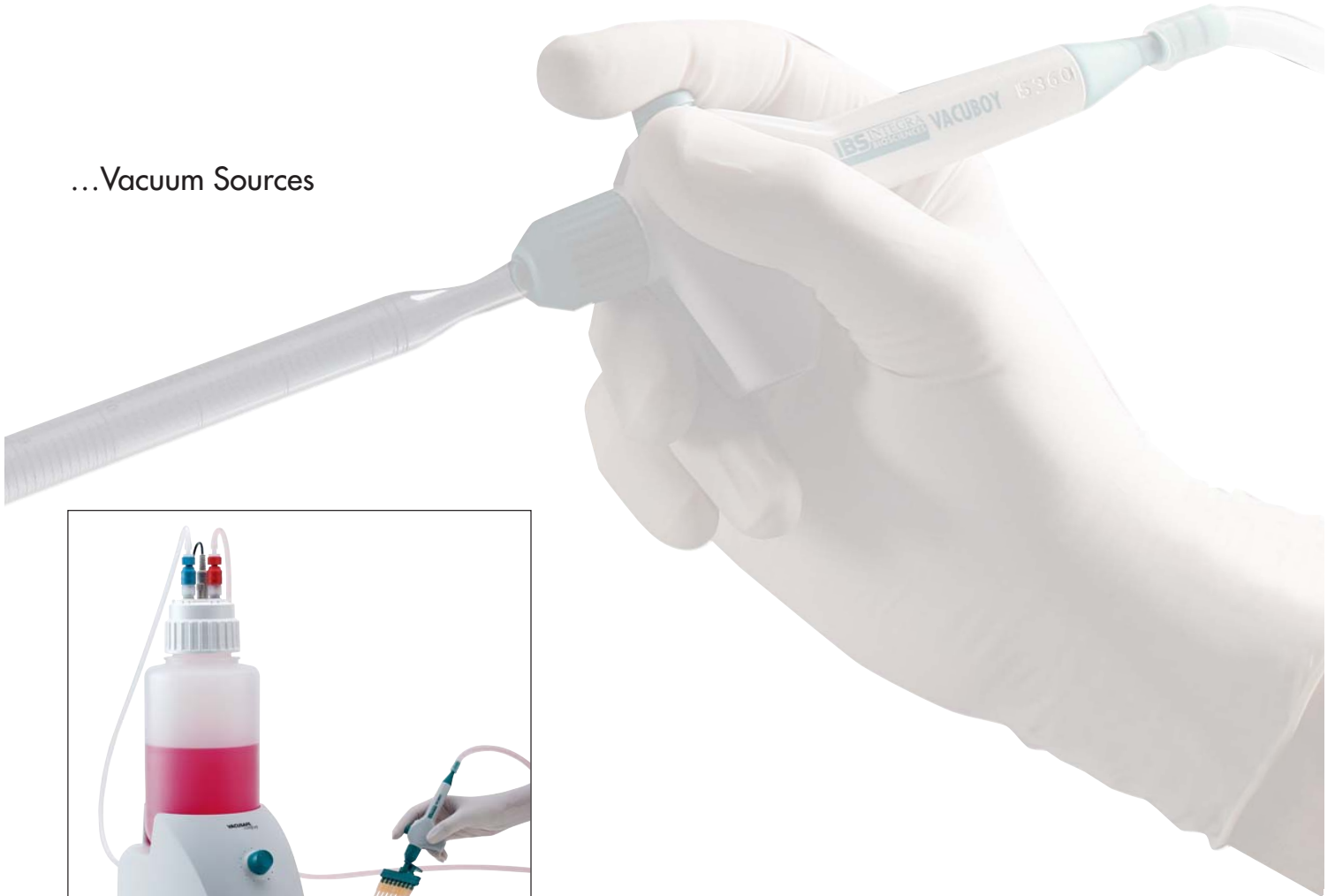


**VACUSAFE** *comfort*




Stand with bottle 4 litre Polypropylene

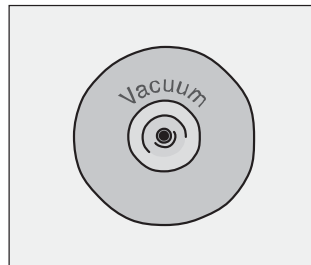

House  
Vacuum System

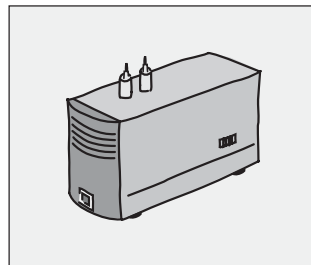

Vacuum Pump

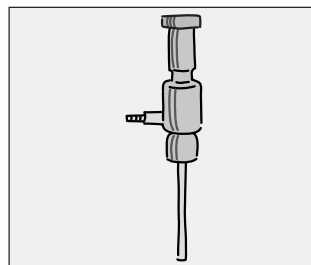

Water Jet Pump

## Ordering Information

Description	Part No.
-------------	----------

**VACUBOY Hand Operator** 155 510

H VACUBOY hand operator

K Rubber adapter for Pasteur pipettes

B 1-channel stainless steel aspiration adapter 40 mm

**VACUBOY Hand Operator Set** 155 500

H VACUBOY hand operator

L Stand for hand operator

K Rubber adapter for Pasteur pipettes

B 1-channel stainless steel aspiration adapter 40 mm

E 8-channel stainless steel aspiration adapter 40 mm

A 1-channel plastic aspiration adapter for disposable tips

I 1-channel plastic aspiration adapter for disposable tips with ejector

G 8-channel plastic aspiration adapter for disposable tips with ejector

## Accessories / Spare Parts

A 1-channel plastic aspiration adapter for disposable tips 155 504

I 1-channel plastic aspiration adapter for disposable tips with ejector 155 526

C 1-channel stainless steel aspiration adapter 150 mm 155 522

D 1-channel stainless steel aspiration adapter 280 mm 155 525

B 1-channel stainless steel aspiration adapter 40 mm 155 502

F 4-channel stainless steel aspiration adapter 40 mm 155 524

G 8-channel plastic aspiration adapter for disposable tips with ejector 155 520

E 8-Channel stainless steel aspiration adapter 40 mm 155 503

Bottle 4 litre Polypropylene 158 370

Lid for Polypropylene bottle with barbed tube fittings 158 372

Hydrophobic membrane filter, non-sterile, Ø 60 mm, 0.45µm 158 015

K Rubber adapter for Pasteur pipettes 155 505

Stand for bottle 4 litre Polypropylene 158 445



L Stand for VACUBOY hand operator 155 501

Tubing 1 800 mm (bottle lid – hand operator) 158 332

Tubing 700 mm (instrument – bottle) 158 331

Tubing set with filter (part no. 158 331, 158 332 and 158 015) 158 342

## VACUBOY Hand Operator Set 155 500

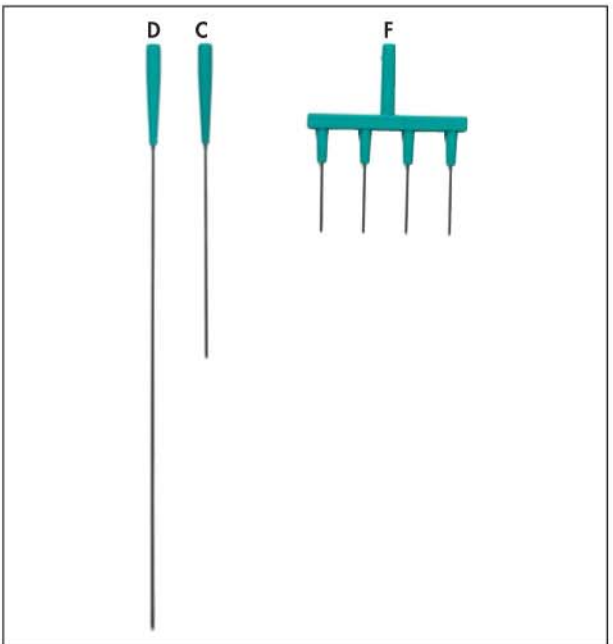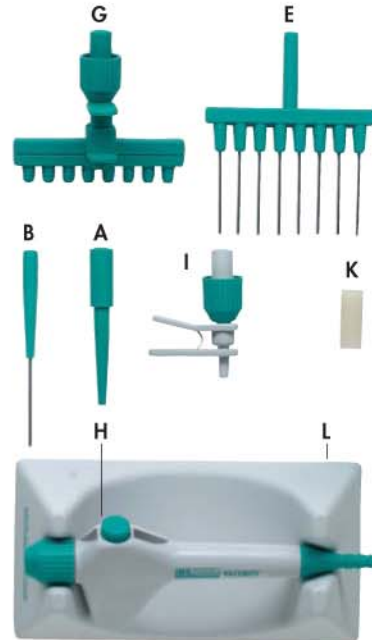

# HYPOXYGEN

v. Define your environment

5111 Pegasus Court, Suite H  
 Frederick, MD 21704  
 877-497-6994  
 www.HYPOXYGEN.us